

**TEXAS PROOF OF CONCEPT AWARD FULL APPLICATION**

# HOME-BASED NON-INVASIVE LUNG COMPLIANCE MONITORING DEVICE

---

Enable non-invasive measuring of lung compliance to monitor and characterize diseases that alter lungs' mechanical properties.

**EDWARD CASTILLO**

Associate Professor, Department of Biomedical Engineering, Cockrell School of Engineering, UT Austin

**HSU-TING KUO**

PhD Candidate, Department of Biomedical Engineering, Cockrell School of Engineering, UT Austin

- Over 16 million U.S. adults have lung diseases affecting the pulmonary parenchyma, costing over \$52 billion annually<sup>1,2,3</sup>.
- Lung compliance (LC)** measures lung's ability to expand at a given pressure, serving as a key indicator of disorders that alter lung's mechanical properties.


	Common pulmonary parenchymal disorders		
	COPD (Chronic Obstructive Pulmonary Disease) <sup>2, 4</sup>	ILD (Interstitial Lung Disease) <sup>3,5</sup>	RP (Radiation pneumonitis) <sup>6</sup>
Prevalence/ Incidence	16M (diagnosed) 24M (undiagnosed)	200K	15-40% of thoracic radiotherapy patients
Cost (USD annually)	50B	2B	10M
5-year survival rate (%)	40-70%	56% (Overall) 41% (IPF) <sup>5</sup>	-
ER visits (in the US annually)	1.5M	13.7%-19.4% of cases	-

- Current lung compliance measuring tools are either **invasive** or **involve radiation exposure from CT imaging**.
- Clinical Unmet Need: No **non-invasive, home-based** solution for regular lung disease monitoring is currently available on the market.

Current LC measuring tools:	Non-invasive?	Precise, consistent pressure control	Frequent monitoring	One-time setup cost	Cost per exam
Mechanical ventilation <sup>1, 2</sup>		✓		\$2.5K – 5K	\$1,522
Pleural manometry <sup>2, 3</sup>				N/A	\$1,038-1,189 (Thoracentesis)
Esophageal manometry <sup>4</sup>				N/A	\$924-1,391
4DCT imaging <sup>5,6</sup>	✓			\$300K -3M	\$6,773
Our Solution	✓	✓	✓	\$700*	\$ 0.75


\* USD 700 covers the cost of the device and the mouthpiece

- Solution: A **non-invasive** system that calculates lung compliance (LC), an important metric in overall lung health, in an **easy-to-use, at-home** package that can be set up and used in minutes.
- Current developmental statue: **Prototyping**
  - Prototype **physical design**
  - Fine-tune **algorithms** for pressure regulation and measurement

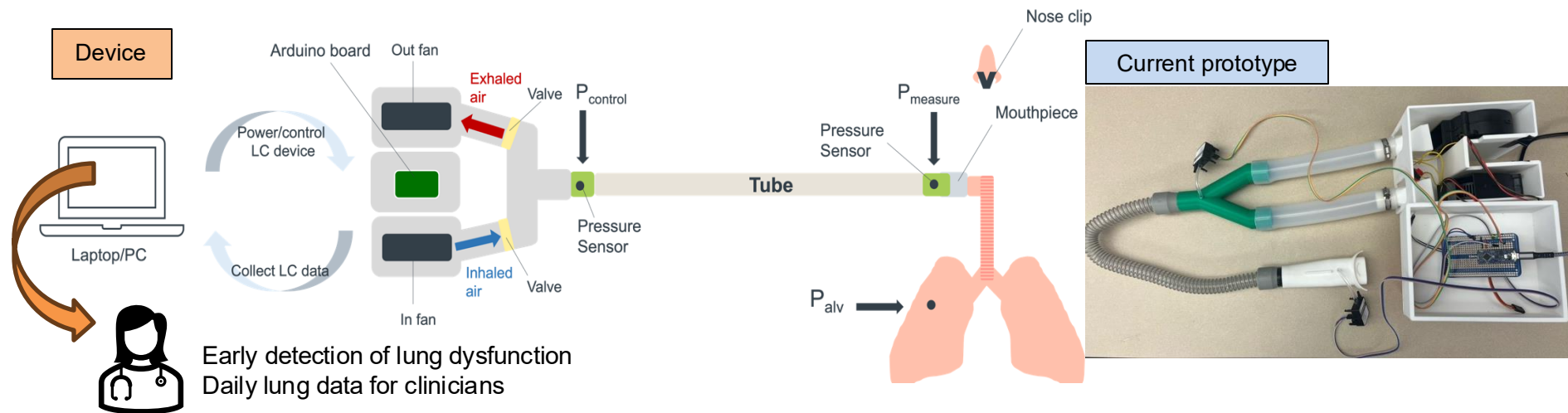

- $LC = \frac{\Delta V}{\Delta P} = \frac{\text{Lung volume change}}{\text{Pressure change}}$  ← Calculated using our numerical pipeline  
 ←  $P_{control}$  across different time points
- Numerical pipeline to estimate airflow rate at all points in the breath cycle
- PID control to regulate pressure during normal breathing cycles

## Numerical pipeline

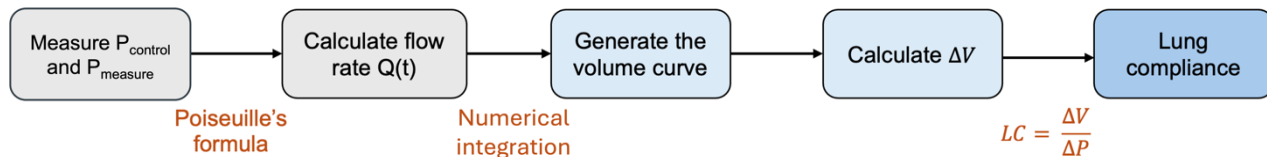

Poiseuille's formula:


$$Q(t) = \frac{P_{control}(t) - P_{measure}(t)}{Res}$$

Numerical integration:

$$Q(t) \cong \frac{dV(t)}{dt} = \frac{V(t + \Delta t) - V(t - \Delta t)}{2\Delta t}$$

Lung compliance calculation:

$$LC = \frac{\Delta V}{\Delta P} = \frac{V(t_1) - V(t_2)}{P_{control}(t_1) - P_{control}(t_2)}$$

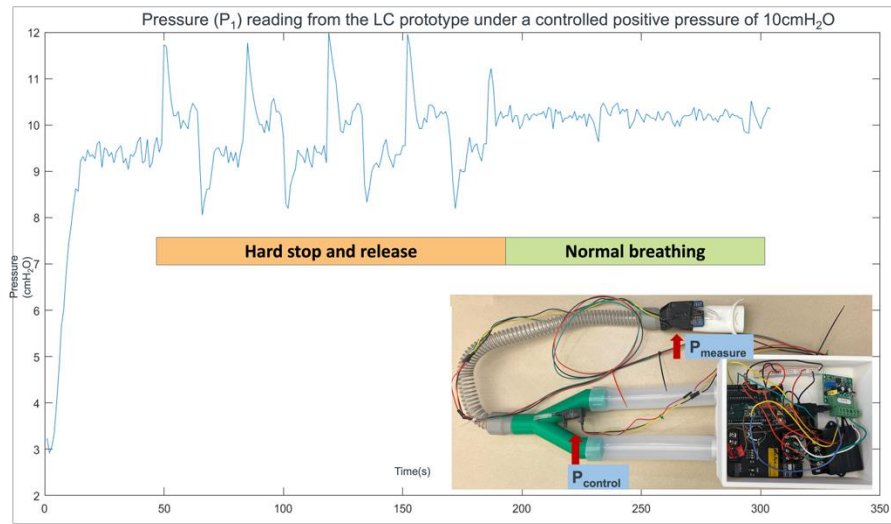

- “Non-Invasive Respiratory Monitoring Market size was over USD 9.09 billion in 2024 and is projected to reach USD 19.87 billion by 2037, witnessing around 6.2% CAGR during the forecast period.” [LINK](#)
- The entry points for our device will be the “Home care” and “Spirometer” segments of the overall market.

## \$350B

**Total Available Market**

People with pulmonary conditions (COPD, ILD, asthma, RP, etc.)

## \$28.7B

**Serviceable Available Market**

The U.S. market

## \$23.25B

**Target Market**

Patients with access to WIFI and healthcare

## 332K


**Users**

Adoption rate = 1%

- The **increasing prevalence of respiratory diseases**, including COPD, asthma, and ILD, is a significant driver for market demand<sup>1</sup>.
  - Global COPD prevalence is projected to reach 600 million cases by 2050, marking a 23% increase from 2020<sup>2</sup>.
- Sanjay Dogra, Director of Respiratory Care at Corewell Health William Beaumont University Hospital, said *“An at-home, user-friendly device for monitoring lung compliance would be valuable for both clinicians and patients...This technology has broad potential applications across a range of lung diseases affecting millions of people worldwide.”*

- Our team has demonstrated significant difference in 4DCT-derived lung compliance (LC) in **healthy** and **IPF** patients<sup>1</sup>.
- We observed **significant decrease** in 4DCT-derived lung compliance over 6 months.

IPF patient

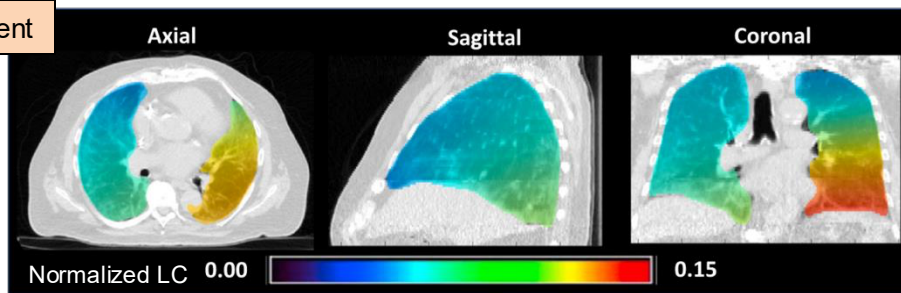

Non-IPF patient

